

# Validator<sup>®</sup> AVS Software 2.1

ACCURACY AND EFFICIENCY WITH SUPERIOR REPORTING

Date Time	Rab1	Rab2	Rab3	Rab4	Rab5	Rab6
24-Feb-2025 12:55:00	21.31	21.83	21.44	20.55	20.48	21.11
24-Feb-2025 13:55:00	21.32	21.83	21.46	20.55	20.48	21.11
	21.32	21.82	21.45	20.55	20.48	21.13
	21.31	21.83	21.44	20.55	20.48	21.11
	21.32	21.83	21.46	20.55	20.48	21.11
	21.32	21.82	21.45	20.55	20.48	21.13
	21.31	21.83	21.44	20.55	20.48	21.11
	21.83	21.46	20.55	20.48	21.11	
	21.82	21.45	20.55	20.48	21.13	
	21.83	21.44	20.55	20.48	21.11	
	21.83	21.46	20.55	20.48	21.11	
	21.82	21.45	20.55	20.48	21.13	

**NEW SOFTWARE FEATURE:  
LIVE HISTORY REPORT**

In today's rapidly evolving digital landscape, where pharmaceutical and biotechnology companies rely on accurate and reliable data, streamlined and flexible reporting is essential.

In addition to the features introduced in AVS 2.0 – such as improved user management and the Common Reporting Tool – the latest version, Validator AVS 2.1, proudly offers enhanced reporting capabilities.

The Historical Live Report provides a real-time table of all collected sensor readings, organized into groups

based on sensors. The collection rate aligns with the sample rate defined in the setup. Users can review data from long-term studies and generate a preliminary Historical Report (PDF or Excel spreadsheet) while the study is still ongoing.

This feature enables users to evaluate data in real time, allowing for early assessment of results. It also offers the flexibility to discontinue long-term studies if early data suggests that the desired outcomes are unlikely to be achieved.


## FEATURES:

- Preliminary History Report (PDF format)
- Preliminary History Report (Excel format)
- Support for KAYE LTR-200
- Active Directory\* support
- Customized user types
- Generated audit and graph reports viewable in the Reports Tile
- Improved data synchronization
- Performance improvements


- Integrated Common Reporting Tool 1.5
  - Merge up to three studies
  - Truncation of large studies to enable graph reports
  - Traceability between qualification studies and associated sensor calibrations
  - Pass/Fail Report
  - Spreadsheet Report
  - Audit Trail Report
  - Parallel report generation with one click




Preliminary History View

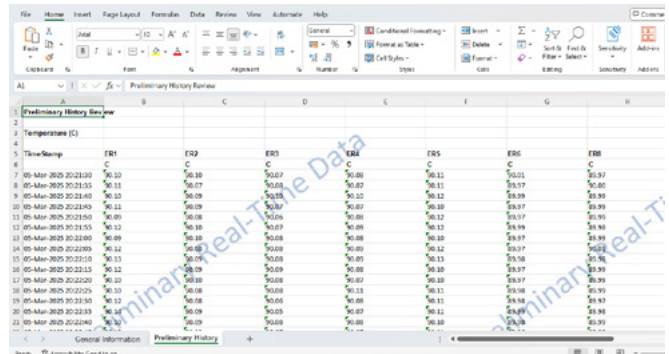

Preliminary History Report (Excel format)

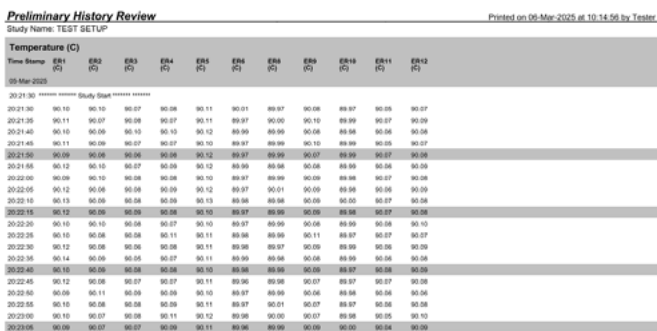

Preliminary History Report (PDF format)

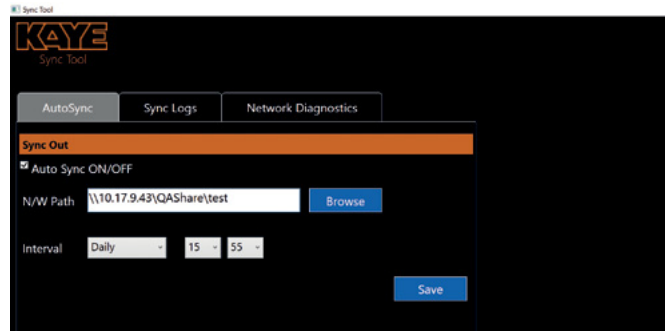

Sync Tool consistently running as a background service

## Kaye Representative contact:

**EUROPE, MIDDLE EAST, AFRICA AND ASIA**  
 Amphenol Advanced Sensors Germany GmbH  
 Sinsheimer Strasse 6  
 D-75179 Pforzheim  
**T:** +49 (0) 7231-14 335 0  
**F:** +49 (0) 7231-14335 29  
**Email:** kaye@amphenol-sensors.com

**USA/AMERICAS**  
 Amphenol Thermometrics, Inc.  
 967 Windfall Road  
 St. Marys, PA 15857  
**T:** +1(814) 834-9140  
**F:** +1(814) 781-7969  
**Email:** kaye-us@amphenol-sensors.com

## Request a demo:

**INDIA**  
 Amphenol Interconnect India Pvt Ltd.  
 Plot no. 6, Survey No.64 | Software Units layout  
 MAHAVEER TECHNO PARK  
 Hitech City, Madhapur | Hyderabad,  
 Telangana – 500081 | T: +91 40 33147100  
**Email:** kaye-india@amphenol-sensors.com

**CHINA**  
 Amphenol (Changzhou) Connector Systems  
 Co., Ltd, Building 10, Jintong Industrial  
 Park, No. 8 Xihu Road, Wujin High-Tech  
 Development Zone, Changzhou,  
 Jiangsu 213164 | T: 0086-519-83055197  
**Email:** kaye-china@amphenol-sensors.com



**Warranty and disclaimer:** The information mentioned on documents are based on our current tests, knowledge and experience. Because of the effect of possible influences in an application of the product, they do not exempt the user from their own tests, checks and trials. A guarantee of certain properties or a guarantee for the proper suitability of the product for a specific, especially permanent application can not be derived from our data. Liability is therefore excluded to that extent permitted by law. Any proprietary rights of third parties as well as existing laws and regulations must be observed by the recipient of the product on his own responsibility.

© 2025 Amphenol Corporation. All Rights Reserved. Specifications are subject to change without notice. Other company names and product names used in this document are the registered trademarks or trademarks of their respective owners.